


The human being is our measure. From the very beginning!
Comfort approved baby scales



M118000 / -01 

Dedicated to health **ADE**

Certified weighing comfort for premature infants, babies and toddlers – also with length measuring system.

We listen and look very closely when it comes to the needs of the smallest patients and the users of our scales. This allows us to pass on the practical experiences of doctors and nursing staff, but also of parents, to our product designers and engineers.

The result combines well thought-out functions with a stable, safe and child-friendly design: Precise, reliable and user-friendly infant and toddler scales for use in surgeries and maternity wards, pediatric clinics and intensive care units.



M118000



M118000 /-01



Scale as film: Please click here!

In focus: **Tested in practice**

**ZERO
TARA**

AUTOMATIC ZERO SETTING

Without any further intervention by the user, the scale automatically finds the zero point every time it is switched on and is immediately ready for weighing.

TARA

Practical it automatically subtracts additional weights and items such as a nappy or a doily from the weighing result.



BATTERY LEVEL ICON

The charge level indication ensures timely battery replacement and thus reliable weighing.

HOLD

REWEIGHING FUNCTION

Second or check weighing at the touch of a button, without having to lift the child off the scale.

AUTO-HOLD FUNCTION

The weighting result is held in the display while the practitioner attends to the patient and can then note the value in peace.

**ON
OFF**


AUTOMATIC SWITCH OFF

The automatic switch-off function, will help to save energy all by itself.

M118000 / -01 

Dedicated to health

ADE

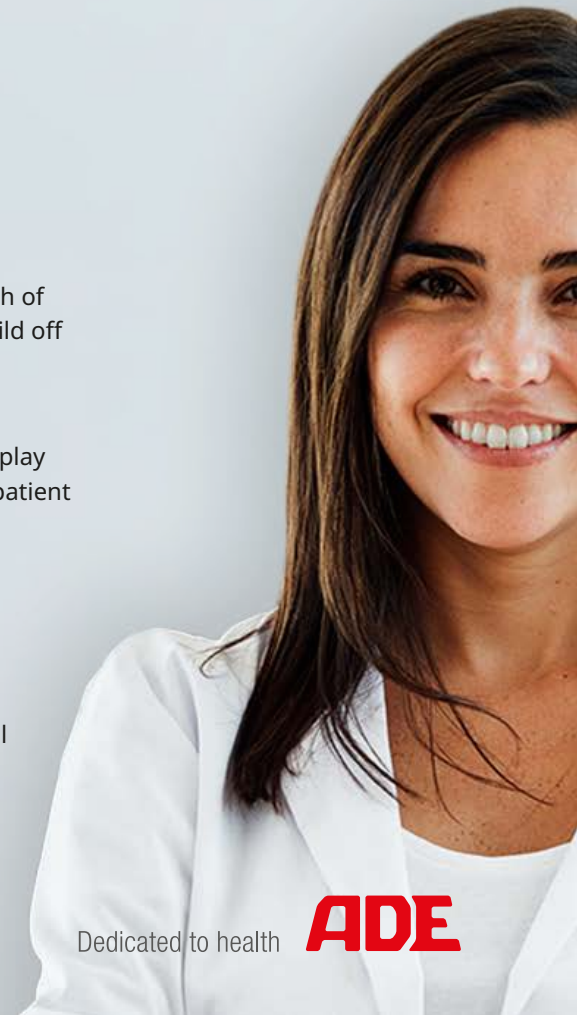

Sophisticated functions



DUAL-RANGE SCALE


Up to a maximum load of 20 kg, these models determine precise results – up to 10 kg with an extra fine division of 5 g, over 10 kg with a division of 10 g.

ECONOMICAL & MOBILE

The scales offer their full range of functions not only in mains operation, but also in mobile battery operation. The 4 (AA) batteries power the scales for up to 4,000 weighings incl. switching on and off.

USER SETTINGS

Extensive selection of individual start settings such as auto switch-off times, auto hold functions on or off, signal sounds and other.

M118000 / -01 


Dedicated to health **ADE**

The technical details

- Capacity 20 kg
- Division 5 g < 10 kg > 10 g
- Backlit display
- Automatic zero setting
- Tara function
- Hold and auto-hold function
- Low battery indication
- Automatic switch off
- Battery operation (4 x AA batteries incl.)
- Mains operation (110-230V) optional
- User settings program
- Internal adjustment program for quick recalibration

Model M118000-01

- Integrated length measurement

M118000 / -01 



The display

- Backlit LCD display for reading in a dimmed environment
- Large high-contrast display
- Robust membrane keypad, protection class IP20
- Intuitive operation via soft-touch keys



M118000 / -01 

Dedicated to health **ADE**

ADE M118000



LARGE WEIGHING SURFACE
with gently shaped tray for a
secure feeling




ROUNDED FORM
without sharp corners
or sharp edges



LIGHT OWN WEIGHT
of 3,2 kg

ADE M118000 

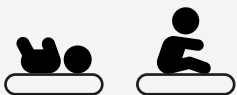


**WEIGHING BABIES
OR TODDLERS** in
laying down or
setting position

DISTURBANCE-FREE WEIGHING
even when holding on to the
edge of the weighing tray

GAP-FREE WEIGHING
surface for easy cleaning



Dedicated to health **ADE**

ADE M118000-01 
with integrated length
measuring function

**FLEXIBLE
POSITIONING**
on the weighing
surface



TWO MEASURING ARMS
for simultaneous
length measuring

Dedicated to health **ADE**

ADE M118000-01 

INDEPENDANT ALIGNMENT
of both measuring arms



EASY SWITCHING
between weighing result and
length measurement



Measuring range from 400 to 800 mm length

Small people. Big task.

Their needs and safety are paramount in the treatment and care of babies. In order for babies to relax in unfamiliar situations, such as being weighed and measured, they need a trustful environment.

That is why we apply great importance to high-quality materials, a pleasant design, stability and quick and easy handling.



M118000 / -01 

Dedicated to health **ADE**

Our quality promise

We offer our customers quality at the highest level. That is why ADE Germany GmbH operates an integrated quality management system which has been certified by TÜV NORD CERT GmbH in accordance with the following standards:

- DIN EN ISO 9001 quality management
- DIN EN ISO 13485 quality management system
for medical devices



Within the EU, weighing instruments for medical purpose and those for determining mass in the practice of medicine must comply with Medical Device Directive 93/42/EEC and Directive 2014/31/EU relating to the making available on the market of non-automatic weighing instruments.



Consulting and contact

+ 49 40 432 776-0
medical@ade-germany.de

ADE Germany GmbH
Neuer Höltigbaum 15
D-22143 Hamburg – Germany
ade-germany.de

M118000 / -01 

Dedicated to health **ADE**