


Column and floor scales

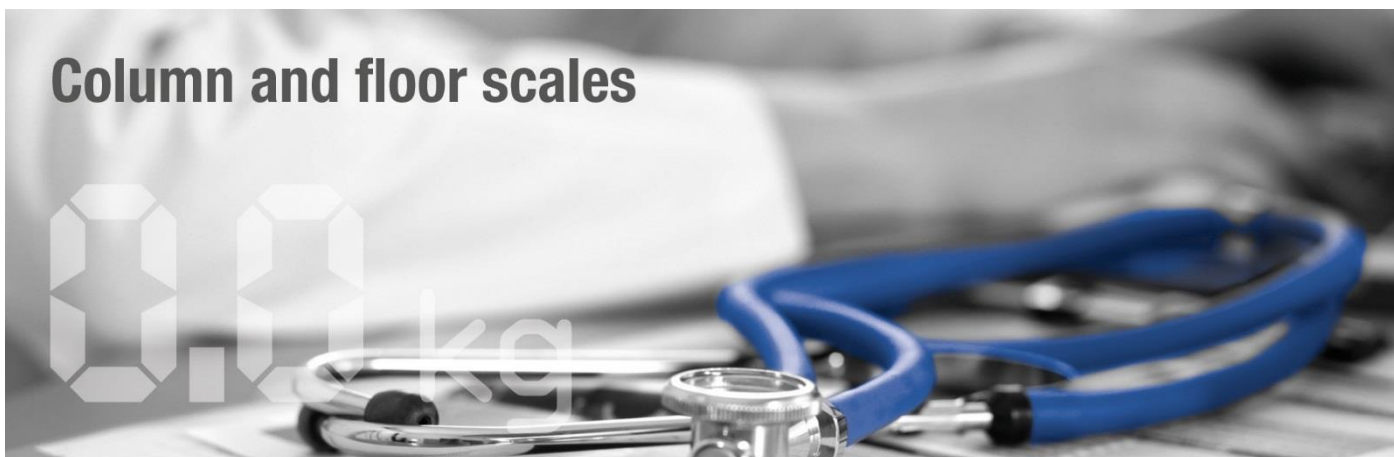

The perfect scale for daily use at Hospitals, Clinics and GPs offices. Strong and sturdy construction with a total capacity of 250 kg and a fine graduation of 100 g. The fully automatic mother and child weighing function will allow determining the weight of new born and small children fast and effortless in the arms of a nurse, parent or field worker just by starting the scale via the tap on function. The display can be turned at a touch of a button towards a helper or the user standing on the scale. Mains (optionally) and battery operated. Accessory ADE MZ10062 transport bag.

Technical specification:

Capacity: 250 kg / 550 lbs
Graduation: 100 g
Dimensions: 400 x 350 x 44 mm
Weight: 3,8 kg

Function: Mother and child weighing, turning display, automatic zero setting, auto switch off, kg/lbs switch. Levelling feet, tap on start function, main all off switch for transport Power supply by 4 x AA batteries life time approx. 4000 cycles - mains optionally.